

# Immedia Positioning Wedge

## SUPPORT AND POSITIONING CUSHIONS WITH NON-SLIP MATERIAL

Immedia Positioning Wedge supports the user on the side after turning, when resting or when performing hygiene or wound care, saving the caregiver from maintaining a static holding position. The Positioning Wedge can be combined with 4WayGlide System, SatinSheet System, MultiGlide or other slide sheets. The cover is a non-slip Polyurethane material that helps to keep the Positioning Wedge in place independent of surface material.



### INTENDED USE

- Support static positioning
- Resting
- Care and personal hygiene
- Can be used alone or in combination with other positioning aids




*The Immedia Positioning Wedge have unequally angled sides, turn the wedge around and use the angle that fits best for the situation.*

# Immedia MultiGlide Glove

## MULTIFUNCTIONAL AID FOR SMALL TRANSFERS

Immedia MultiGlide Glove can be used in many different manual handling situations. The glove makes it easier to place the hands between the user and mattress or cushion. Small re-positioning adjustments, movements and smoothing out the bed sheets are common activities for which the gloves are very useful,

and will provide an improved solution for the user and the caregiver. MultiGlide Glove can also be used for skin tension release.

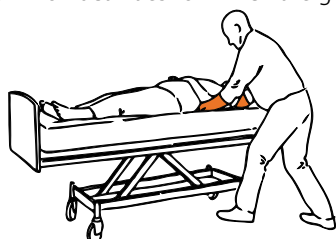

### INTENDED USE

- Positioning
- Small re-positioning adjustments
- Hand check of anti-decubitus cushions and mattresses



*Immedia Padded MultiGlide Glove.*

*Immedia MultiGlide Glove in waterproof nylon.*

### DESCRIPTION

### WIDTH X LENGTH

### ITEM NO.

### SPECIFICATIONS

#### Immedia Positioning Wedge

	CM	INCH	
~ Small	30 x 20 x 53	11.8 x 7.8 x 20.9	IM53NS
~ Large	30 x 20 x 80	11.8 x 7.8 x 31.5	IM53/80NS
~ Cover for Positioning Wedge, Small	30 x 20 x 53	11.8 x 7.8 x 20.9	IM5300NS
~ Cover for Positioning Wedge, Large	30 x 20 x 80	11.8 x 7.8 x 31.5	IM5300/80NS

Can be wiped with disinfectant.



#### Immedia MultiGlide Glove

	CM	INCH	
~ MultiGlide Glove waterproof nylon (1 pair)	20 x 50	7.9 x 19.7	IM35
~ MultiGlide Glove Padded (1 pc)	15 x 35	5.9 x 13.8	IM8000701

#### WASHING INSTRUCTIONS

Do not use fabric softener or bleach.



#### Immedia Positioning Wedge

MATERIAL: Soft foam  
COVER: Nylon with PU-coating, zipper



#### MultiGlide Glove, padded

INSIDE: Nylon

OUTSIDE: Nylon

#### MultiGlide Glove, waterproof nylon

MATERIAL: Nylon with polyurethane coating.

